

Septihol™ WFI Sterile Alcohol Wipes combine knitted polyester wipes with a ready to use (RTU) solution of 70% (v/v) USP Isopropanol and 30% USP Water For Injection (WFI) to provide a presaturated wipe. The right-sized packaging system offers a convenient way to introduce low-particulate wipes to your production environment, while reducing product waste.

Septihol WFI Sterile Alcohol Wipes are part of the Septihol family of products. Each product is designed to simplify your validation requirements while providing multiple packaging options.



Septihol™ WFI Sterile Alcohol Wipes

Ready-To-Use Alcohol Wipes

Features

70% (v/v) USP Isopropanol and 30% USP Water For Injection alcohol solution

Specifically formulated for easy evaporation

Right-sized packaging

Knitted polyester wipe with sealed edges

Double bagged and gamma irradiated, with every lot USP sterility tested

Each lot tested for endotoxin levels of <0.25 EU/mL, meeting the Water for Injection (WFI) specification

Benefits

Industry best practice for cleaning and rinsing

Ideal for use on critical surfaces

Reduced waste of wipes

Low particulate wipes to achieve compliance with current Annex 1 contamination control strategies


Facilitates introduction into critical areas

Suitable for use in critical areas

Applications

Surface cleaning of:

- Countertops
- BSL hoods
- Carts and other equipment

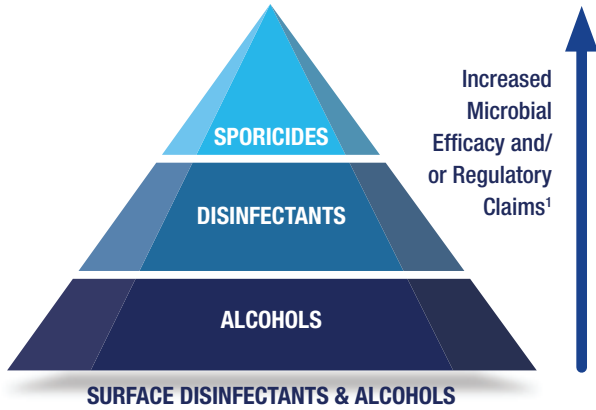

¹Products that fall into the categories at the bottom of the pyramid are most frequently used and are generally not sporicidal. Progression up the pyramid indicates stronger performance overall and a broader spectrum of claims.

When following cleanroom cleaning procedures and aseptic techniques, it is necessary for personnel who perform critical activities to clean work surfaces with the proper agent to prevent contamination. It's important to use sterile products that evaporate quickly and completely, and do not leave a residue. The application methods of cleaners used in any microbial control program are critical and should be carefully selected. They must be designed to balance validation needs for consistency and repeatability with operational needs for application and performance. Septihol WFI Sterile Alcohol Wipes contain the optimal concentration of 70% (v/v) USP Isopropanol for cleaning work surfaces, offering a simple, repeatable application method for areas where mopping or spraying are difficult. Our Technical Service group can help customize your contamination control program to make full use of this product.

Septihol WFI Sterile Alcohol Wipes Ordering Information

Description

Pouch of fifteen (15) pre-saturated 12 by 12 inch (30.5 x 30.5 cm) wipes with 427 mL of chemistry
12 pouches per case
(total of 180 wipes per case)

Part Number (United States and Canada)

624915