

## Custom Washing Rack Scanning Services

Maximize Washing Efficiency with Precision Scanning

### Tailored Racks, Trusted Performance

Improve efficiency, reduce lead times and streamline operations with precise washer rack models. STERIS's on-site scanning team offers a unique service, collaborating with your facility to design custom washing racks fitted to your specific parts and equipment, ensuring accurate and effective washing cycles.



# How it Works

Our comprehensive approach starts with a thorough on-site 3-D data collection of your parts and equipment. This leads to a determination of appropriate rack options, detailed proposals, and the design of your washing rack. We validate every aspect through rigorous testing and compliant documentation, ensuring a tailored result that meets your unique washing needs.



## 1. 3-D Data Collection

The STERIS scanning service team visits your facility to collect and record measurements.



## 2. Evaluation

After collecting the data, we determine appropriate options, adjusting your rack design to your specific needs.



## 3. Proposal

Following the evaluation process, our team will provide you a detailed proposal and pricing for the custom washing rack. We can also conduct a Computational Fluid Dynamics (CFD) flow dispersion analysis upon request or if our engineering team deems it necessary.



## 4. Design


Once you approve the proposal, the order will be sent to our design and manufacturing facility, where our team will design your rack. We will then send you a package for approval before production begins.



## 5. Validation

After you approve the final rack design and it goes into production, we will provide you with GMP-compliant documentation and offer you the option to supply your parts for thorough coverage testing and cleaning validation using your specific soil, if provided.

STERIS's scanning services team handles every aspect of data collection and rack design at your facility. We precisely measure each desired piece of your equipment and design custom racks that securely accommodate and support them, enhancing the efficiency of your operations.

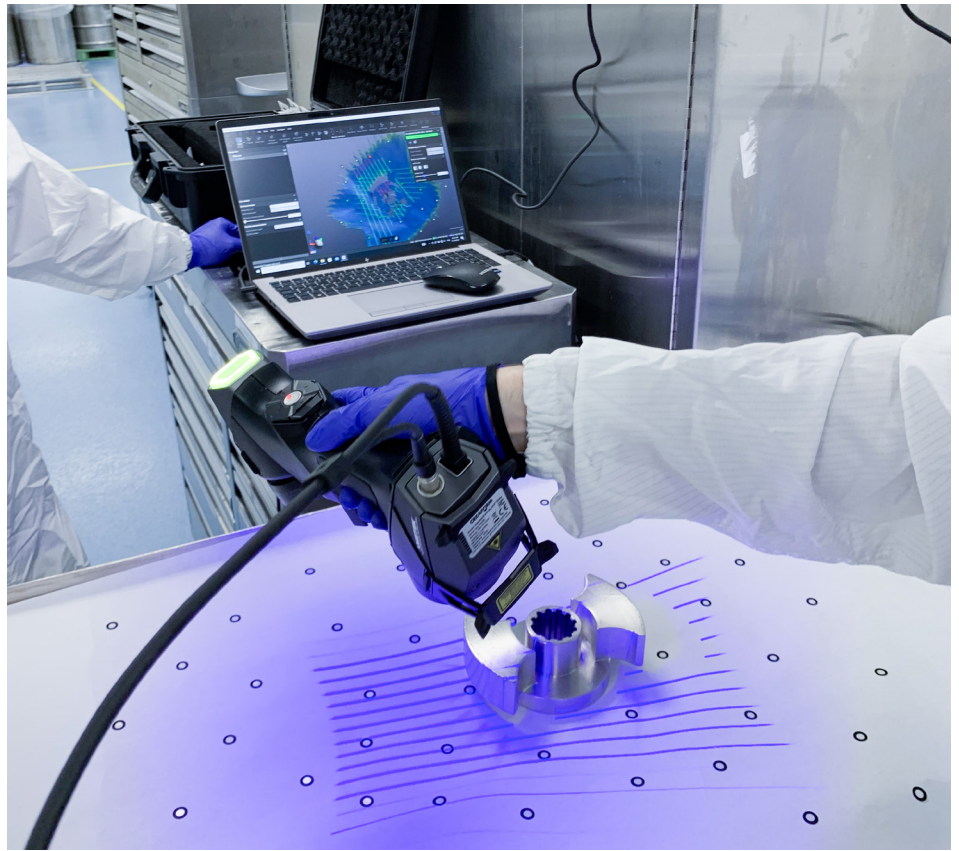

Benefit	Service Feature
<p><b>Improve Efficiency</b></p>	<p>A custom rack ensures that all your parts are correctly loaded and washed each time. This not only reduces lead times but also enhances your facility's operations by optimizing the cleaning process.</p>
<p><b>Save Time</b></p>	<p>This service streamlines the data collection of unique parts and equipment, eliminating the need for you to gather this information yourself.</p>
<p><b>Ensure Accuracy</b></p>	<p>Our scanning technology provides accurate measurements of every part, optimizing washing cycles, improving wash performance and extending equipment lifespan.</p>
<p><b>Lower Costs</b></p>	<p>This approach to custom washing racks ensures thorough coverage for your parts and equipment, enhancing cleaning effectiveness while minimizing discrepancies and operational costs in your washing cycles.</p>
<p><b>Streamline Validation</b></p>	<p>Our custom washing racks come with comprehensive documentation to support compliance efforts.</p>

